


BeneSalus Technologies Smart Pillbox

Our Smart Pillbox is a Bluetooth Low Energy (BLE) Daily Pillbox for reminders for 4 different times through the day.

Along with the MedMonitor app, it forms a powerful tool that helps you take charge of your health with better compliance.

In the following pages we present some of the features of the box

BeneSalus Technologies Smart Pillbox

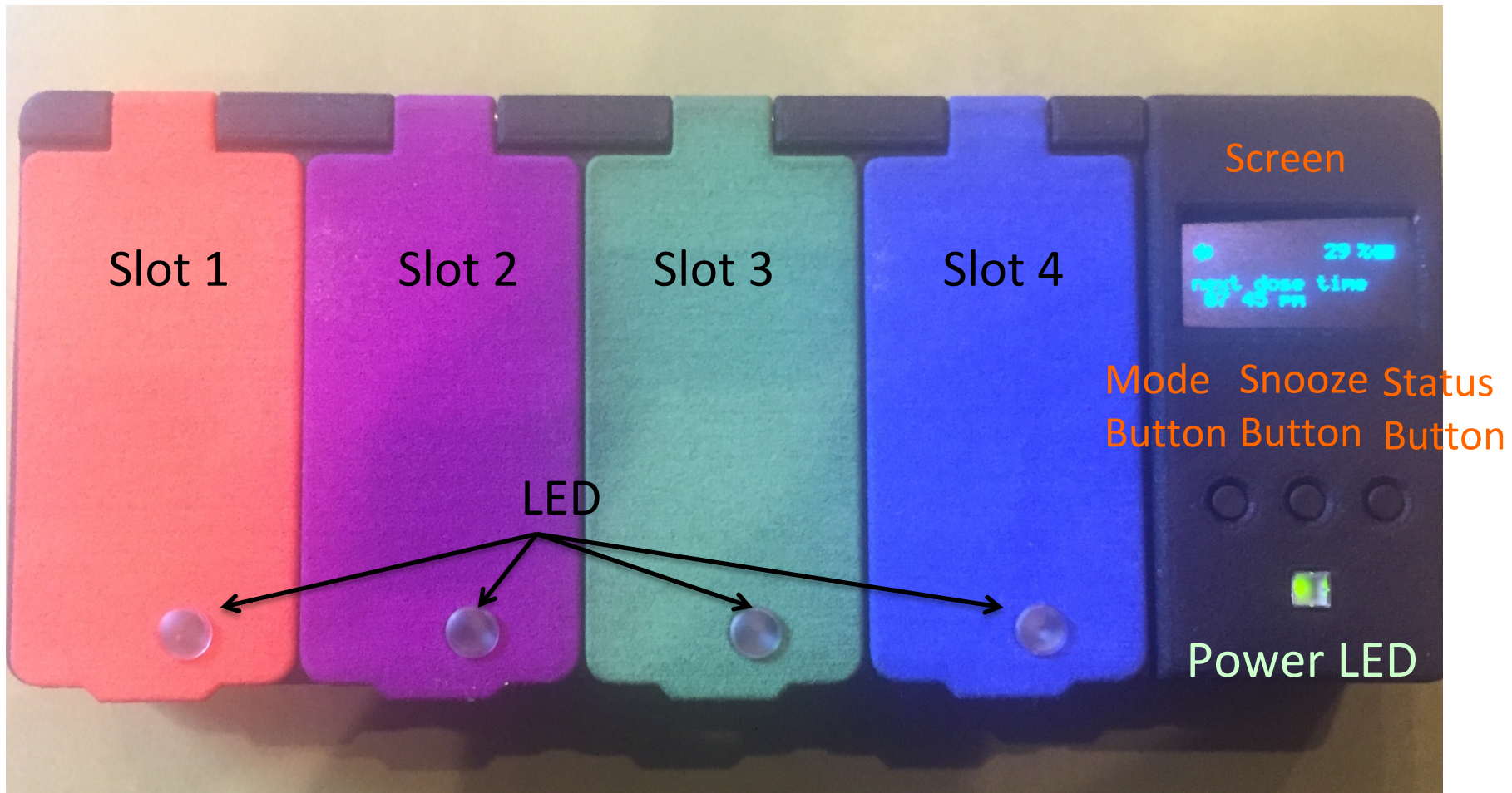

Smart Pillbox Features

- Pillbox connects with MedMonitor App via Bluetooth Low Energy (BLE)
- The App available on Apple App and Google Play stores, configures the Pillbox with for Medications and Dose times information.
- The Pillbox will remind you when to take your medications when it is time for a particular dose
- Pillbox features both an audible Buzzer and a Vibrator
- Pillbox can be charged via USB or a Wireless charger. A full charge should last for a week of use

Smart Pillbox in Action

- The box reminds you when it is time to take your medications
- It does it as follows
 - It will sound its buzzer
 - The LED for the correct pill slot with medications to be taken will blink
 - The Screen will display the name of the medication(s) to be
- All you have to do is open the Indicated slot and take your medication. *Please remember to keep the lid open for at least 5 seconds for the box to register the event*
- The pillbox registers the event and updates the MedMonitor App and your Circle of Care with the time the medications were taken when it next syncs with the App

Loading the Pillbox

When done with all doses for the day, you are ready to reload the box with medications for the following day. Do that as follows

- Start by Pairing the Pillbox with the App.. Bring the Pillbox close to your Smartphone and open the app. The box will pair automatically
- Proceed to the Settings page of the app. The Pillbox controls on the App Settings page will remain hidden till the Pillbox is connected to the App. You'll have to wait a few seconds for that to happen.
- Put the Pillbox in **Load Mode** using the Load Mode toggle control
 - LEDs for all slots that need to be loaded will blink. Note that if your app is configured for 2 doses then only the first two slots will blink
- Open the lids for all the slots that need to be loaded
- Take care to put all medicines in their corresponding slots. The slot numbers are labeled on the lid
- When done loading, close all the slots
- Remove the Pillbox from Load Mode. The LEDs will stop blinking and the box is ready for the next day

Pillbox Buttons

There are three buttons on the pillbox as show in Figure 1. There functions are explained below

- The three buttons are Mode, Snooze (Snz) and Status (Stat)
- The Mode button allows the box to be in their an Audible buzzer mode or a Silent Vibrate Mode
- Snooze allows you to snooze the Dose Alarm for 5minutes
- Status Button Identifies the box as your and gives an update on missed doses and when the next does time is

Pillbox Buttons

- To power the Pillbox ON, press and hold the Mode button till the display turns on.
- The Display will show “Connect to Phone”
- Make sure Bluetooth is enabled on your phone
- Open the MedMonitor App on your phone and the Pillbox will automatically connect to the App
- For more info and help visit

www.benesalustech.com/support.php

Pairing the Pillbox with MedMonitor App

Congratulations on getting the Smart Pillbox

To power the Pillbox ON, press and hold the Mode button till the display turns on.

Mode button is the right most button under the screen

The Display will show “Connect to Phone”

Make sure Bluetooth is enabled on your phone

Open the MedMonitor App on your phone and the Pillbox will automatically connect to the App

For more info and help visit www.benesalustech.com/support.php

support@benesalustech.com